

Flower Orthopedics Shoulder Solutions




The Flower Proximal Humerus Plate Construct consists of a high-strength stabilization plate designed to best accommodate the anatomy and pathology of the proximal humerus. The anatomic curvature of the Flower Proximal Humerus plate is designed for efficient placement and can be bilaterally applied to right and left fractures. The Flower Proximal Humerus Plate accepts variable angle locking, non-locking and compression screws in diameters 2.7mm to 4.0mm.

Features & Benefits


- Anatomic contour facilitates reduction and limits need for bending
- Proximal suture holes aid in reduction and tuberosity fixation
- Plates accept multiple screw options for patient and pathology variability
- Screws sit below the the surface of the plate, reducing risk of soft tissue irritation
- Variable angle screw placement up to 30° cone

Proximal Humerus Plates

- Address complex fractures, fracture dislocations, osteotomies and nonunions of the proximal humerus
- Anatomic bilateral design available in 5 lengths
- Screw holes are positioned for optimal subchondral support and fixation of marginal fragments

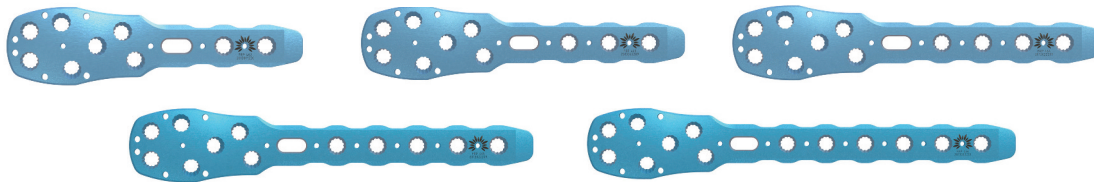

Screw Options	Length
3.5mm Variable Angle Locking Screw	10-60mm
3.5mm Variable Angle Non-Locking Screw	10-60mm
3.5mm Variable Angle Compression Screw	10-40mm
4.0mm Variable Angle Compression Screw	10-60mm

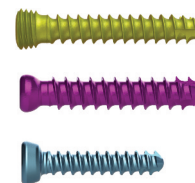

INDICATION FOR USE: The Flower Small and Medium Implants set is intended for use for internal fixation of fractures and reconstruction of bones, including the scapula, olecranon, humerus, radius, ulna, pelvis, distal tibia, fibula, hand and foot in adults and for long bones in adolescents (12-21) in whom the growth plates have fused. Examples of these internal fixations and reconstructions include compression fractures, intra-articular and extra-articular fractures, displaced fractures, osteotomies, non-unions and mal unions. This system can be used for palmar, ventral, dorsal or orthogonal application.

©2015 Flower Orthopedics. All rights reserved. FlowerCube, FlowerCarriage, IndicationCube, ApplicationCube, IndicationCarriage, E-Kit, FlowerCart, Ready-for-Surgery, FlowerPatch, FlowerFlo and FlowerPot are trademarks of Flower Orthopedics. Please refer to package insert for description, indications, contraindications, warnings, precautions, and other important information.

Flower Orthopedics Corporation • 100 Witmer Road, Suite 280 • Horsham, PA 19044 • T: 877.778.8587 • F: 215.394.8904 • www.flowerortho.com