

Flex-Thread™ Ulna Intramedullary Nail System

The Flex-Thread Ulna Nail System is an innovative solution for the treatment and fixation of ulna fractures and osteotomies.

SYSTEM FEATURES

Innovative Flex-Thread Design

- Flexible titanium nail shaft
- Atraumatic blunt nose tip design
- Self-contours to native anatomy
- Screw-like insertion assists in maintaining reduction

Proximal Nail Fixation


- Three divergent screw trajectories optimize rotational stability
 - 30° proximal screw hole to engage the olecranon
 - Coronoid specific screw trajectory
 - Transverse distal screw trajectory

Distal Nail Fixation


- Provides secure point-contact fixation
- Multiple thread diameters
- Reverse cutting flutes

ZERO Profile Interlocking Screw Design


- 3.5mm titanium self-tapping
- Minimizes soft tissue irritation
- All screws thread into nail for improved fixation

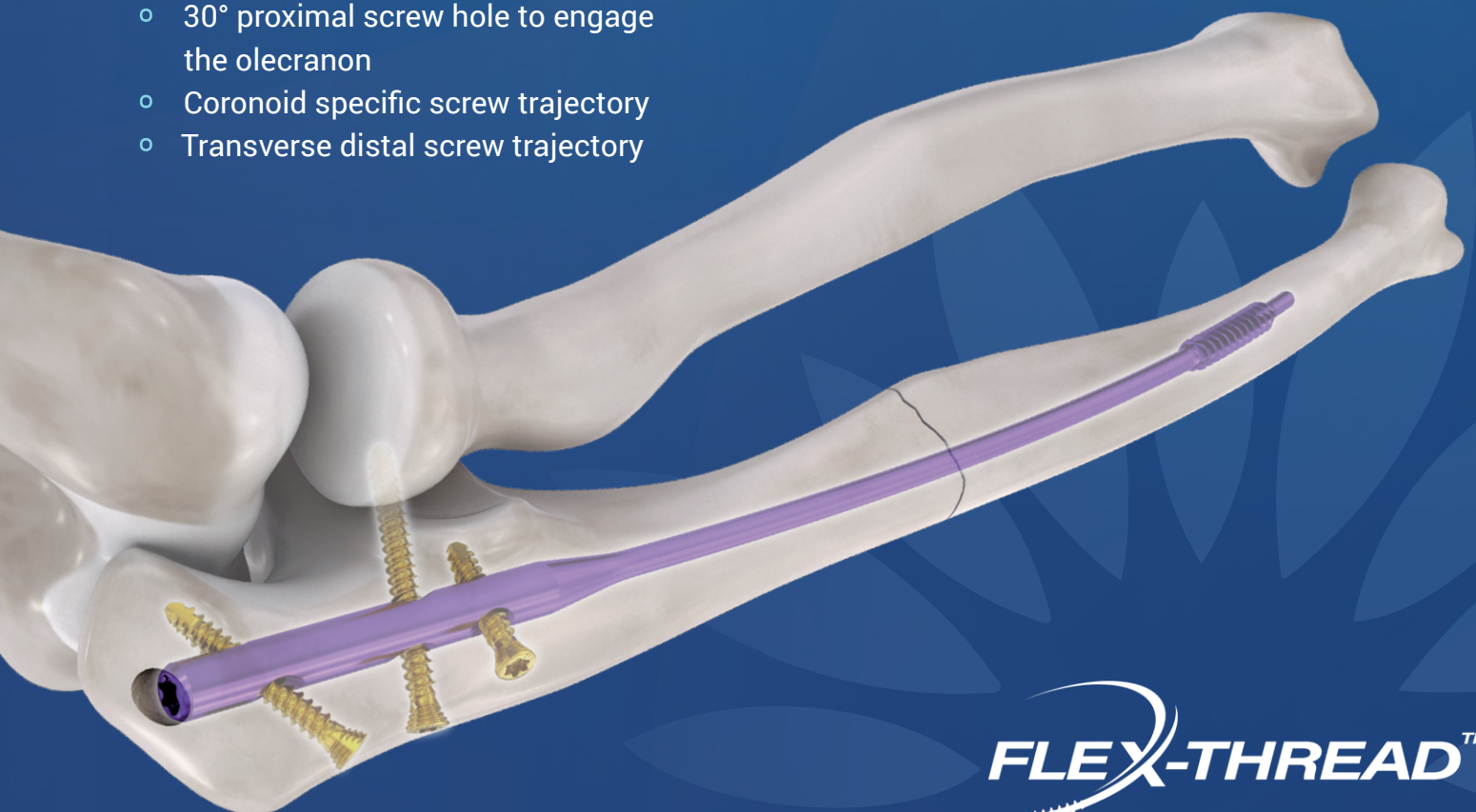

FLEX-THREAD™

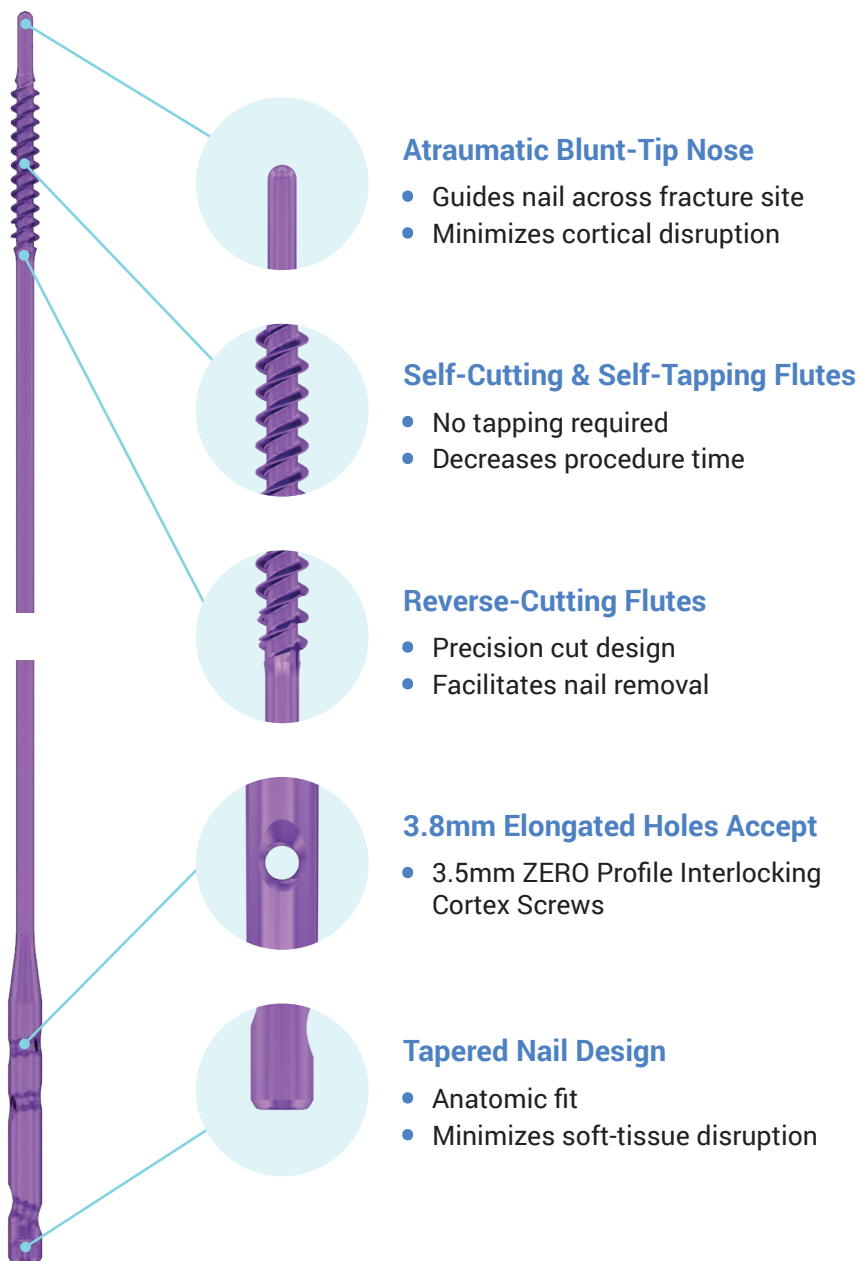

Nail Diameter		Nail Length	
4.5mm	210mm	240mm	270mm
5.5mm	210mm	240mm	270mm
6.5mm	210mm	240mm	270mm

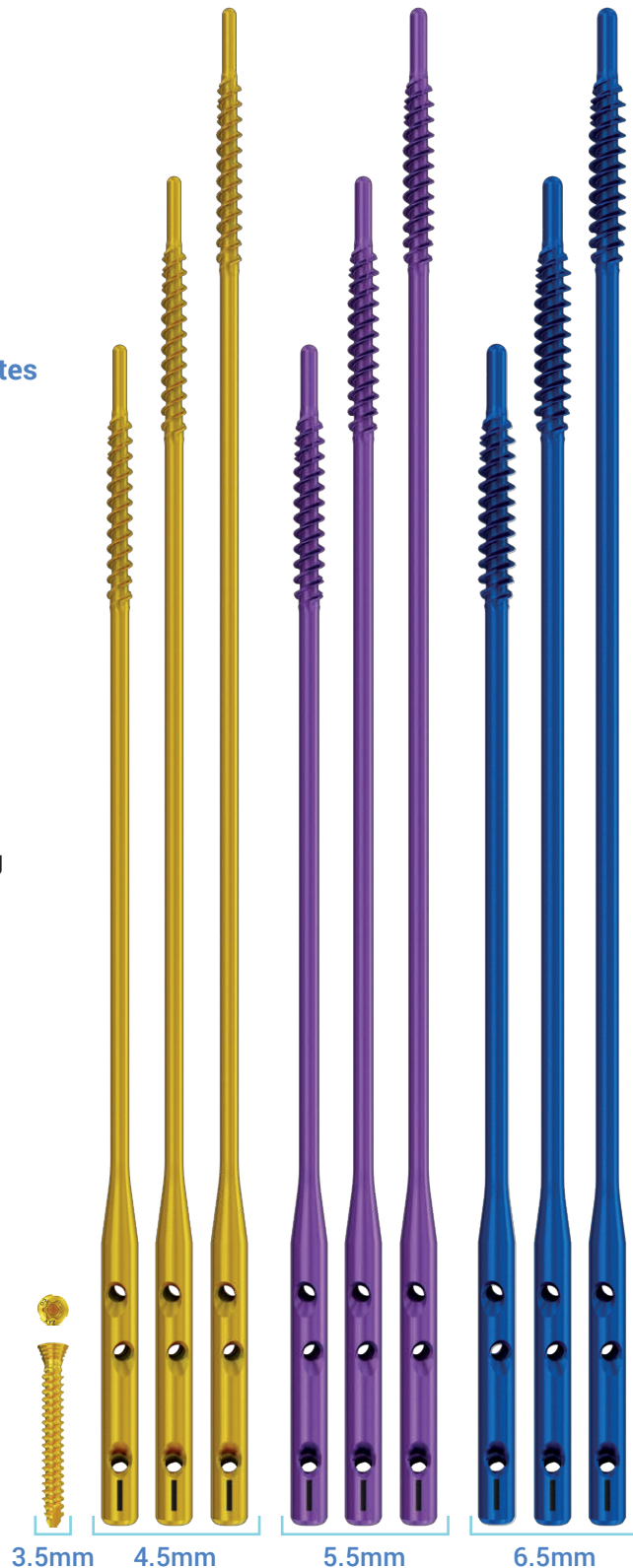

Conventus Flower Orthopedics
100 Witmer Road, Suite 280, Horsham, PA 19044
877.778.8587

[ConventusFlower.com](https://www.ConventusFlower.com)